

ionic Nū Lift. Rejuvenation System

Younger Looking Skin For Life.



The ionic **Nū-Lift™ Rejuvenation System** by Essona Organics is a powerful, lightweight, multi-purpose device with state-of-the-art galvanic micro-current technology designed to help draw out impurities and to increase the absorption potential of the skin for immediate results.

Expensive treatment at your doctor or beauty salon is no longer necessary, and you just need about 8-16 minutes to help draw out impurities and prevent the appearance of wrinkles, improve skin tone, elasticity and vibrancy.



Science Behind the Technology

The skin tissue is composed of epidermis as the outer layer, dermis, and hypodermis that contains fat. Skin cells renew every 28 days. The protective membrane in the epidermis layer can isolate skin from the harmful environment and substances and protect it, but at the same time, it may hinder the skin from absorbing nutritious substances.

The ionic **Nū-Lift™ Rejuvenation System** uses both micro current and galvanic current with specific polarity to help draw out the impurities on a deeper level. Meanwhile, the gentle current stimulates the skin with its thorough massage, which supports the development of healthy skin tissues and strengthens collagen and elastin fibers. Thus skin becomes delicate and soft, fresh and natural, full of elasticity, with a younger-looking glow.

The ionic Nū-Lift™ Skin Rejuvenation System


Essona
organics™


www.essona.com

Product Specifications:

Power Supply	2 PCS AAA Battery	Measurement	135 X 38 X 23 mm
Working Current	10 ~ 50 mA	Weight	43 g
Rated Voltage	3V		FDA Certified
Basic Features	Four step beauty system: Cleansing, Massage, Nutrient and Lifting		
	Micro current & Galvanic technology with both positive ion and negative ions		
	Vibration micro-massage with vibration sensor		
	Auto stop function with 4 minutes for each step		
	Easy to use, portable and can be used anytime and anywhere		

Components & Function

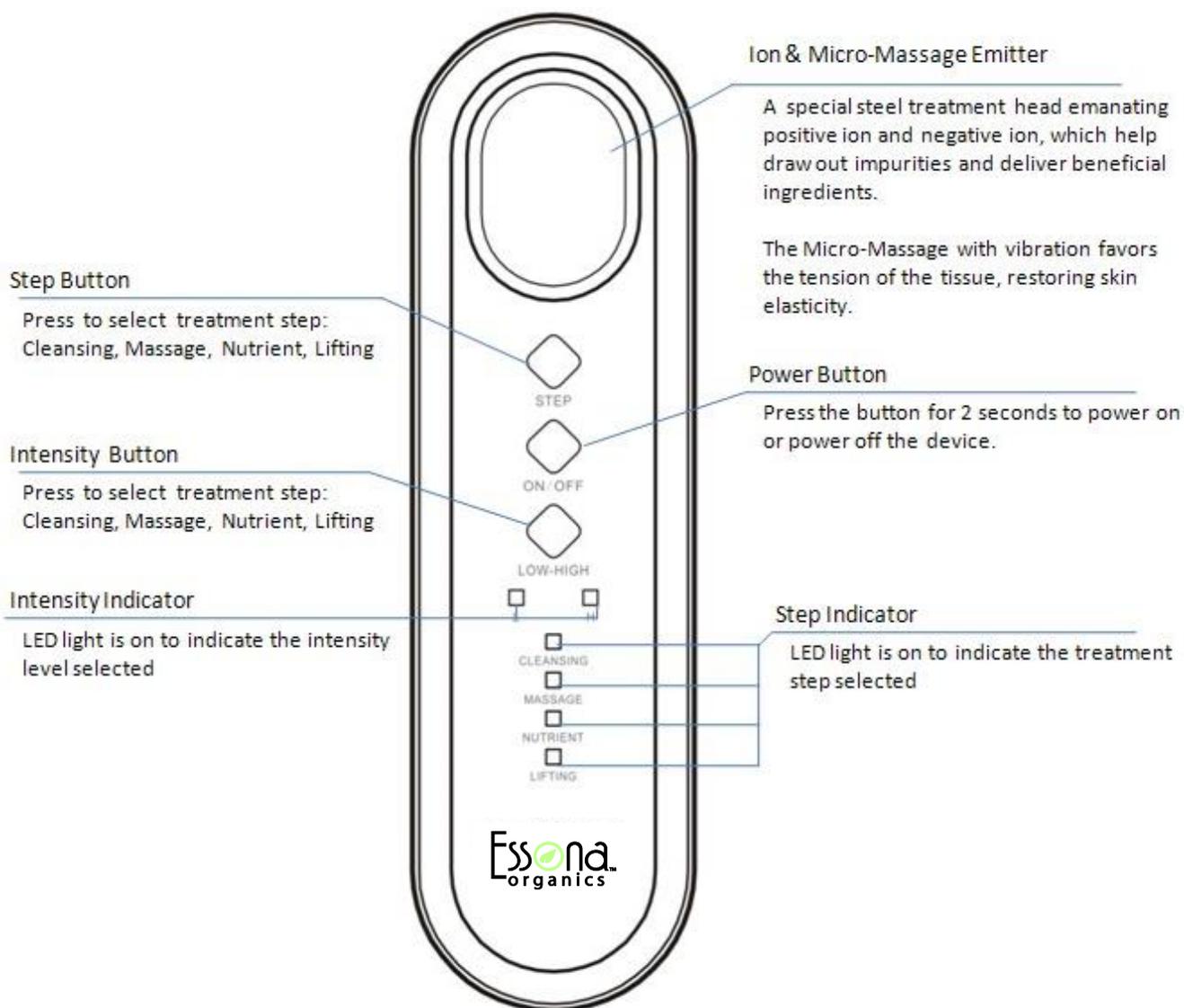

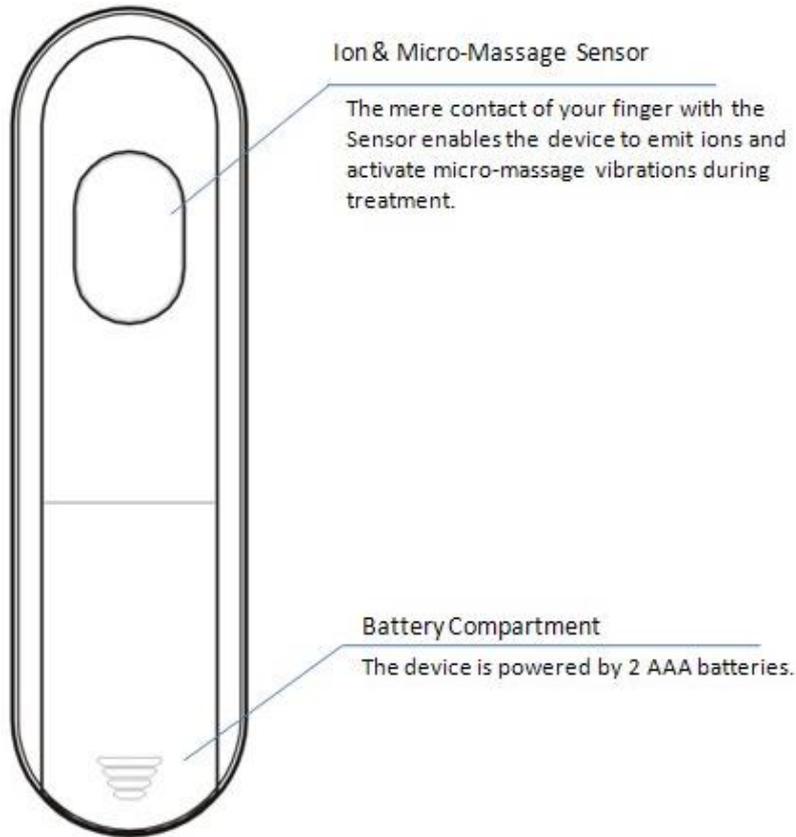

The **ionic Nū-Lift™ Rejuvenation System** is designed with four treatment steps for a professional spa experience at home. Each step should last 2-4 minutes, depending on the desired level of results. After 4 minutes, it will stop automatically.

Uses 2 AAA batteries.

Do not use oily creams while using the unit. Use water-soluble lotion or conducting gel to get the best conductivity and results. Use the Essona Organics skin care treatment regimen that is specifically designed for use with the **ionic Nū-Lift™ Rejuvenation System**.


Step
1

Cleansing

Dirt inside
the pores

Dead cells

Grease &
Sebum

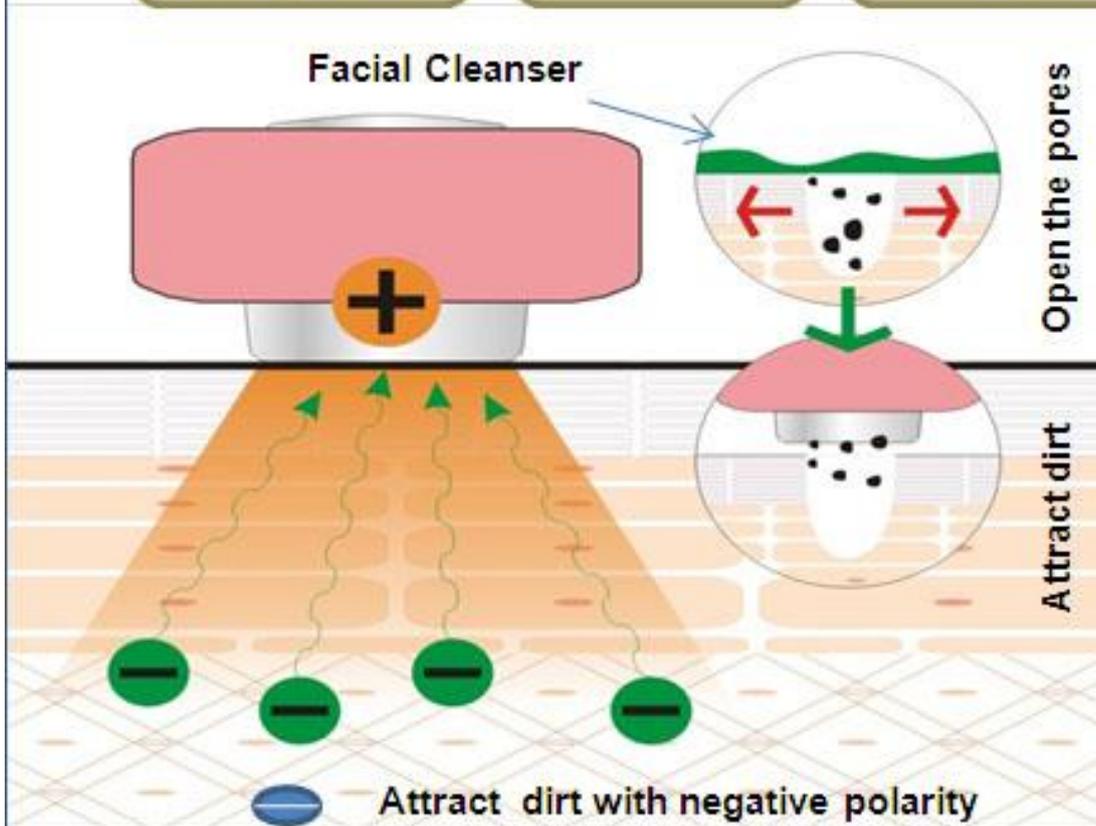

Attract impurities inside the pores, deep cleansing skin!

The positive ion attracts the blackhead and dead cells deeply from the pores, removing excess sebum deposits, dirt, debris and impurities.

Use: Apply facial cleanser, gently massage the face with machine in circular movements from inside to outside, from down to up.

Step
2


Balance

Tired skin

Dull skin

Dry skin

Massage products

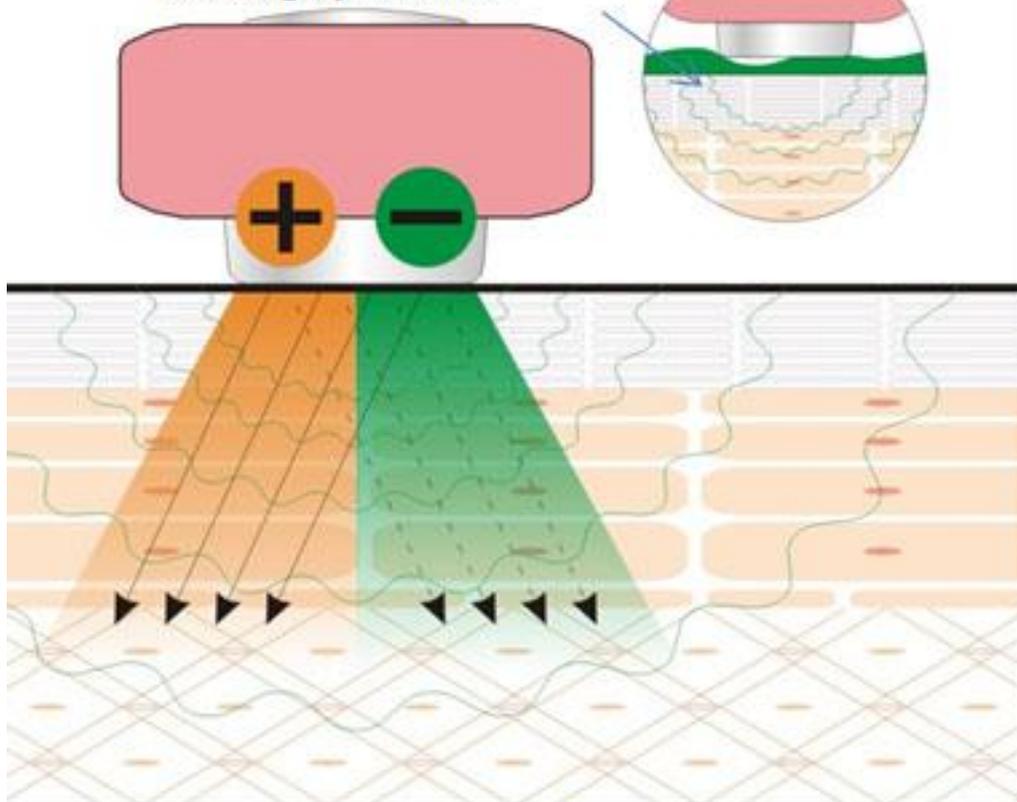

Soothing current stimulation, refresh skin!

The positive and negative ion stimulate the skin alternatively, promoting blood circulation and metabolism of skin cells, making the skin more healthy and glossy.

The Step 2 prepares the skin for Step 3.

Use: Apply massage products, gently massage the face with machine in circular movements from inside to outside, from down to up.


Step
3

Nutrient

Dull skin

Dry skin

Big pores

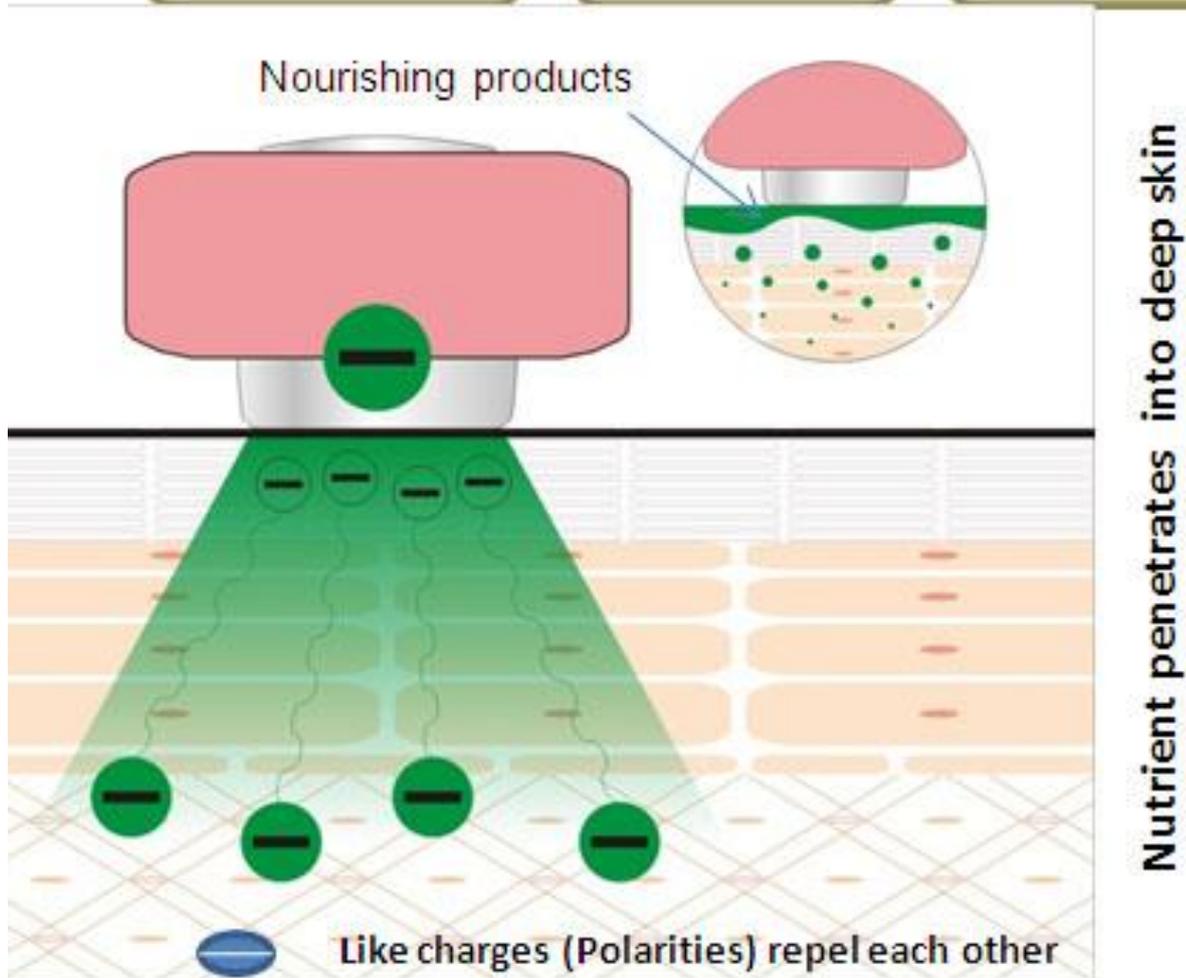

Greatly drive nutrient into deep skin!

The negative ion makes the water-soluble ingredients and nutrients into the skin, moisturize skin and make skin smooth, vivid and tender.

Use: Apply nourishing products, gently massage the face with machine in circular movements from inside to outside, from down to up.


Step
4

Lifting

Sagging skin

Wrinkle skin

Aged skin

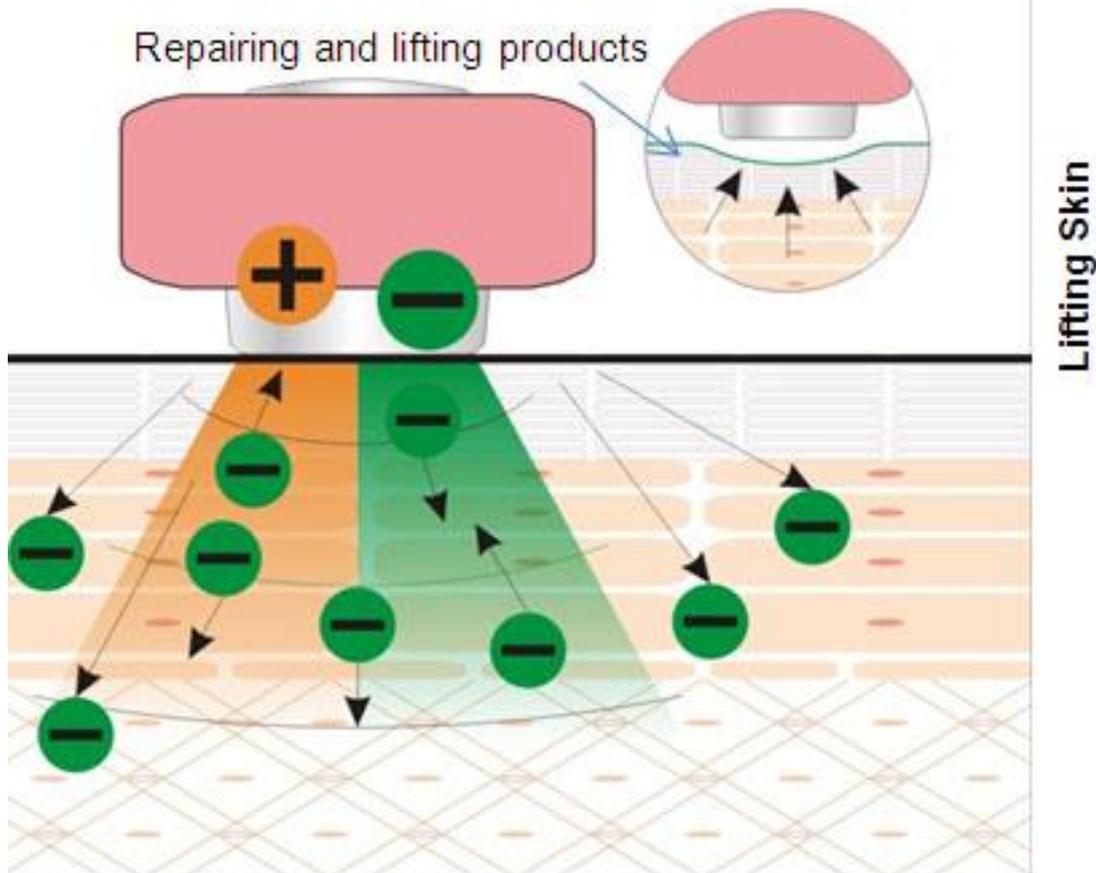

Lifting sagging and dropping skin and reshape facial lines!

The positive and negative ion stimulate skin alternatively, reduces face muscles tension, enhances the connective tissue, erases deep wrinkles and keeps skin tight, lifting.

Use: Apply repairing and lifting products, gently massage the face with machine in circular movements from inside to outside, from down to up.

Four Step System to Reveal Younger-Looking Skin:

Step 1. Cleansing – Draw out impurities inside the pores, deep cleansing skin.

Remove makeup and use Spring Dew Facial Cleanser to wash your face. Rinse well. Apply NuFusion Treatment Gel generously on your whole face or specific area of concern. Switch unit to “Cleaning” setting. Choose intensity level (start with the “low” setting). Gently work the unit in a circular motion from the center of your face outward. Avoid neck and eyelids. The positive ions attract blackheads and dead cells deeply from the pores like a magnet, removing excess deposits of dirt, debris and impurities better than cleansing alone. (2-4 minutes).

Step 2. Balancing - Soothing current stimulation to refresh skin.

Switch unit to “Balance” setting. Choose intensity level (start with the “low” setting). Apply NuFusion Treatment Gel if necessary to keep your skin lubricated. Gently work the unit in a circular motion from the center of your face outward. Avoid neck and eyelids. The positive and negative ions stimulate (massage) the skin alternatively, promoting blood circulation and increased metabolism of skin cells, making the skin more healthy and glossy. (2-4 minutes).

Step 2 prepares the skin for Step 3.

Step 3. Nutrition – Nourish your skin at a deeper level.

Switch unit to “Nourish” setting. Choose intensity level (start with the “low” setting). Apply NuFusion Treatment Gel if necessary to keep your skin lubricated. Gently work the unit in a circular motion from the center of your face outward. Avoid neck and eyelids. The negative ions help to invigorate and improve the penetration potential of the healthy nutrients in the Essona Organics Skin Care products, thereby leaving your skin healthy, smooth and glowing! (2-4 minutes).

Step 4. Lifting- Support sagging and drooping skin and reshape facial lines.

Switch unit to “Lifting” setting. Choose intensity level (start with the “low” setting). Apply NuFusion Treatment Gel if necessary to keep your skin lubricated. Gently work the unit in a

circular motion from the center of your face outward. Avoid neck and eyelids. The positive and negative ion surge stimulates skin alternatively, reduces facial muscle tension, enhances the underlying connective tissue, erases fine lines and wrinkles and keeps skin smooth and tight.

(2-4 minutes).

*After you are done with your session, rinse your face and apply SeaMist Toner and NoLines Face Serum and NoLines Eye Serum. Wait a few minutes and finish with Ancient Secret (lifting), Wine Rose (anti-aging), Pearl Glow (whitening, brightening,) Royall Lotus (nourishing), Seatopia (hydration) Herbal Bliss (healing), Clearal (acne) or your favorite Essona Organics product.

***Note: You do not have to do all 4 steps at once. You can spread the step throughout the week. The 2 most important steps are cleansing and balancing. Make sure that you keep your face well lubricated with NuFusion Gel to increase conductivity and reduce friction. The only product you should use while using the unit is NuFusion Treatment Gel.**

The ionic Nū-Lift™ Rejuvenation System is designed to be used with Essona Organics Skin Care products.



Inner Health, Outer Beauty

www.essona.com